

# DryMax Soft


DryMax Soft is a sterile, conformable, superabsorbent dressing for exuding wounds. Its properties provide an easy-to-use and reliable solution for both patient and clinician. DryMax Soft absorbs and retains exudate in the dressing. The exudate is absorbed vertically by the dressing and retained in the core, forming a gel, thus minimizing the risk of maceration and maintaining a moisture balanced wound environment.

The seamless cover simplifies positioning over the wound. The design allows the dressing to be applied with either side facing the wound. This feature also makes it possible to stack dressings on each other if extra absorption capacity is required.


The dressing remains dry and comfortable for the patient. Simultaneously, it contributes to a moist wound environment. The hydrophilic contact layer transports fluid effectively away from the wound and distributes it evenly into the dressing's core constructed by our DryMax Technology.

The ultra-thin and highly conformable features allow the dressing to be used on difficult external body parts. DryMax Soft is designed to work well together with compression therapy.


## Application Guide

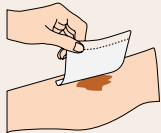

The dressing must be placed close to the wound, or placed in close combination with a suitable dressing.

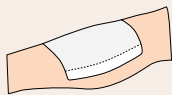

Secure the dressing with suitable bandage, film or tape.

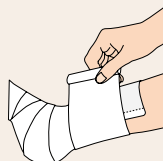

DryMax Soft can be used in conjunction with compression therapy.

Item	Product size cm	REF	Pcs/box	Absorption-capacity*
DryMax Soft	5x5	F60125/10	10	20 ml
DryMax Soft	7.5x7.5	F60011/10	10	60 ml
DryMax Soft	10x10	F60013/10	10	130 ml
DryMax Soft	12.5x12.5	F60123/10	10	220 ml
DryMax Soft	10x15	F60012/10	10	300 ml
DryMax Soft	10x20	F60015/10	10	330 ml
DryMax Soft	15x20	F60014/10	10	400 ml
DryMax Soft	20x20	F60016/10	10	600 ml
DryMax Soft	20x25	F60018/10	10	780 ml
DryMax Soft	20x30	F60017/10	10	960 ml
DryMax Soft	20x40	F60124/10	10	1300 ml
DryMax Soft	37x56	F60019/10	10	3000 ml

The maximum free swell absorption capacity of DryMax Soft is 1.8 ml per cm<sup>2</sup> of the product core.

\*According to international standard EN 13726

\*\*compared with conventional wound dressings



- **Cost effective**
- **User friendly**
- **Increased patient comfort**
- **Improved wound healing\*\***



**DryMax® technology**