

DryMax Easy


DryMax Easy is a sterile, conformable, superabsorbent dressing for exuding wounds. Its properties provide an easy-to-use and reliable solution for both patient and clinician. DryMax Easy absorbs and retains exudate in the dressing. The exudate is absorbed vertically by the dressing and retained in the core, forming a gel, thus minimizing the risk of maceration and maintaining a moist and balanced wound environment.

The soft tack function helps the dressing to adhere to the skin and will work as an extra hand during the dressing application. The dressing can be applied with either side facing the wound. This feature also makes it possible to stack dressings on top of each other if extra absorption capacity is required. The seamless cover simplifies positioning over the wound.


The dressing remains dry and comfortable for the patient. Simultaneously, it contributes to a moist wound environment. The hydrophilic contact layer transports fluid effectively away from the wound and distributes it evenly into the dressing's core constructed by our DryMax Technology.

The ultra-thin and highly conformable features allow the dressing to be used on all body parts. DryMax Easy is designed to work very well together with compression therapy.


Application Guide

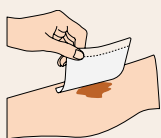

The dressing must be placed close to the wound, or placed in close combination with a suitable dressing.

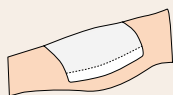

Secure the dressing with suitable bandage, film or tape.

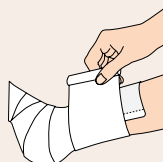

DryMax Easy can be used in conjunction with compression therapy.

Item	Product size cm	REF	Pcs/box	Absorption-capacity*
DryMax Easy	10x10	F60033/10	10	170 ml
DryMax Easy	10x20	F60035/10	10	500 ml
DryMax Easy	20x20	F60036/10	10	1000 ml
DryMax Easy	20x30	F60037/10	10	1500 ml

*According to international standard EN 13726



- **Cost effective**
- **User friendly**
- **Increased patient comfort**
- **Improved wound healing***

*Compared to conventional wound dressings



DryMax® technology